SONIC TRAVEL

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PL

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USER'S MANUAL

EN


Dear Customer, thank you for purchasing **EVOREI SONIC TRAVEL**. The product you have obtained is a sonic toothbrush. This device is used to brush your teeth and maintain oral hygiene. **Before using the device, please read the instructions below below.**

01. SAFETY MEASURES

- If you do not understand the instructions for use or safety measures listed in this manual, do not use the device. Incorrect use of the device may damage it.
- Do not disassemble the product or replace any of its parts. This may result in damage to the product, electric shock and void the warranty.
- Do not use accessories other than those recommended by the manufacturer.
- Do not use flammable or explosive substances to clean the product.
- This product is only used for brushing teeth and massaging gums. Do not use the device for any other purpose.
- Stop using the device if it is clearly damaged.
- Do not place the product in water. This can lead to a malfunction and may result in fire, electric shock or explosion.
- Keep the device away from water containers such as a bath or sink.
- Do not replace the battery of the device when your hands are wet.
- Do not use old batteries in your device.
- Pay attention to the poles of the battery when placing it in the device.
- Do not place the device in rooms with high moisture or high temperature, such as bathrooms.
- Remove the battery from the device if you will not use it for an extended period of time to avoid the risk of fire.
- If the toothpaste used contains peroxides or bicarbonate (often present in tooth whitening pastes), wash the device thoroughly after use to avoid destroying its plastic casing.
- If your gums bleed or you feel discomfort for at least a week of use, consult your dentist before continuing to use the device.
- The device is not intended for use by persons with limited physical, cognitive or mental abilities or who lack experience or knowledge, including children. If such a person wishes to use the device, they must be supervised or instructed by the person responsible for them.

- The device is not a toy – keep it away from children's range.
- The device is used for personal hygiene and should not be used on many people in dental offices.

02. DESCRIPTION OF THE DEVICE

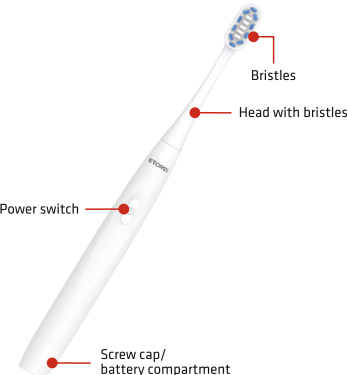03. PACKAGE CONTENT

- Toothbrush
- Head with bristles
- User's Manual

04. POWER

The device needs 1 piece of AAA battery for proper operation. Unscrew the lid at the bottom of the device to access the battery compartment. Place the AAA battery according to the diagram below.

Place the cover back in place and screw it firmly in place. The bump on the lid should be under the bottom of the arrow on the housing before screwing, and under the dot after screwing.

05. USING THE DEVICE

- Put the head with bristles on the metal core sticking out of the device body.
- Squeeze the toothpaste onto the bristles.
- Gently press the toothbrush to your teeth and gums.
- Press the power switch to turn on the brush and start cleaning your teeth. Thoroughly clean every part of your teeth.
- When you're done brushing your teeth, turn off the toothbrush again by pressing the power button.
- After cleaning your teeth, wash the toothbrush bristles under running water. Remove the remaining water from the bristles and use delicate material to wipe the handle of the device.
- When you are not using the brush, keep it upright.

NOTE: To avoid splashing toothpaste foam, turn off the device before removing it from your mouth.

06. MAINTENANCE

- CLEANING THE BRISTLES:**
 - Wash the brush head and its bristles after using the device and let them dry naturally.
 - Do not press or pull the bristles out with your fingers. This can deform the toothbrush and destroy it.
- CLEANING THE DEVICE HANDLE:**
 - Wipe the handle with delicate material and let it dry naturally.
 - Do not submerge the handle under water.
 - Do not use water above 40°C or acidic cleaners to clean the appliance.

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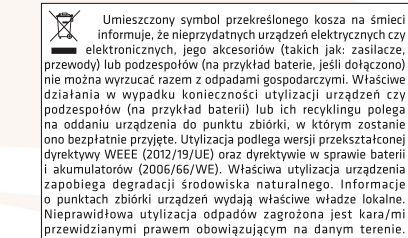
Producent:
KGK Trade sp. z o.o. sp. k.
Os. Uroczce 12
31-953 Kraków
Polska

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
09. WARRANTY/ CLAIMS

This product is covered by a 24-month warranty. The warranty conditions are available on the website: <https://xblitz.pl/warranty/>

Complaints should be submitted with the complaint form which can be found at <http://reklamacje.kgktrade.pl/>

Contact details and the address of the service can be found on the website: www.xblitz.pl

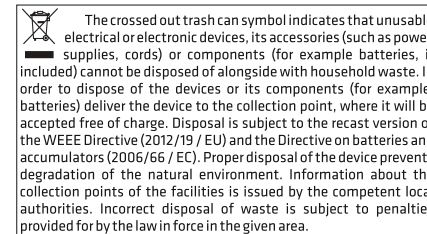Manufacturer:
KGK Trade sp. z o.o. sp. k.
Os. Uroczce 12
31-953 Cracow
Poland

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